


TPN Teaching Diagram



Wash your hands and prepare supplies.

1

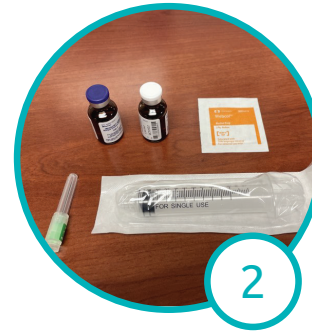

Remove caps from multivitamins, 1 blue top, 1 white top, clean tops with alcohol pads. Withdraw multivitamins from the vials into the syringe after carefully attaching the needle to the syringe.

2



Clean port on TPN bag with an alcohol pad and carefully inject multivitamins into the TPN bag. Dispose of needle and syringe in sharps container. Mix TPN bag by rotating.

3



Spike TPN bag with tubing.

4



Power on the pump and open the latch. Place tubing as directed and close the latch. Do not connect tubing to patient yet.

5

6

Select yes to prime tubing. Once line is filled and consistent drip is coming out, select stop.

7

Clean the clave (connector) leading to the patient and connect the TPN tubing.

8

Ensure all clamps are open and press start.